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
USER MANUAL

LED DÉCOLLETÉ MASK | BEAUTY TOOL WITH
LIGHT THERAPY

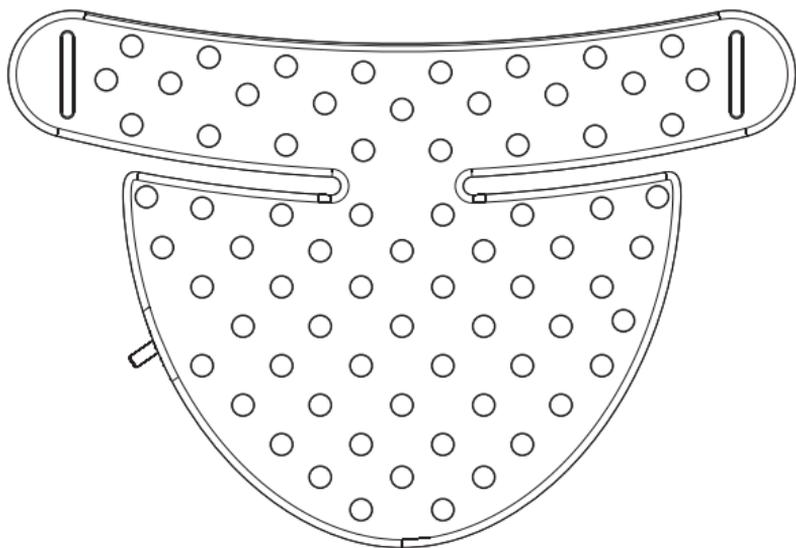FIRMING

REVITALIZING

BALANCING

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USER MANUAL

Enjoy every moment with your LED Décolleté Mask | Advanced Beauty Tool that includes light therapy from Rosental Organics!


Before using your new tool, please read this manual carefully and keep it for future reference.

Our LED Décolleté Mask is a home-use device designed to promote collagen production, improve circulation and skin elasticity, and reduce inflammation, blemishes, and pigmentation.


To ensure correct usage of this product, please read all safety precautions carefully before using. These precautions will help you use this product correctly to avoid personal injury or damage to yourself or others. The safety precautions are listed below. Please comply with the corresponding regulations!

PARTS

Remote Control

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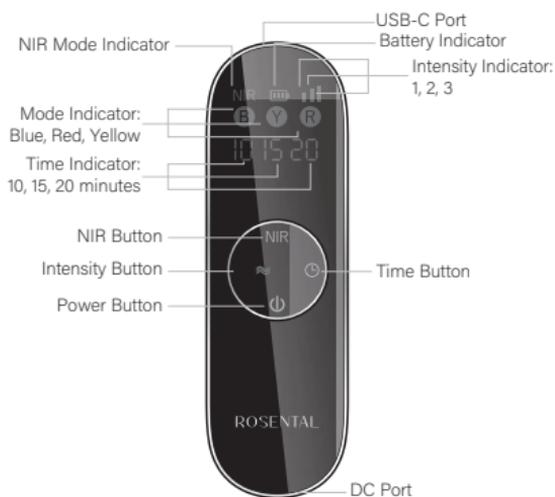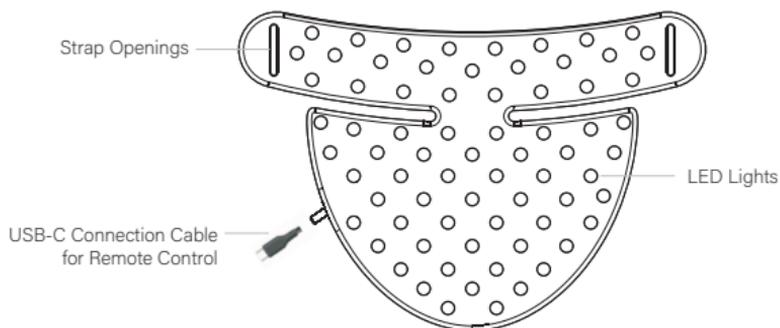
ASSEMBLY

STEP 1

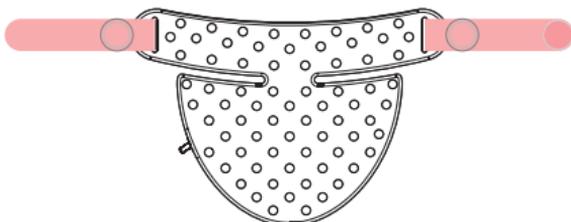
Attach the straps to the mask by threading them through the openings and securing them inside with the Velcro.

STEP 2

Remove the protective film from the mask.

STEP 3

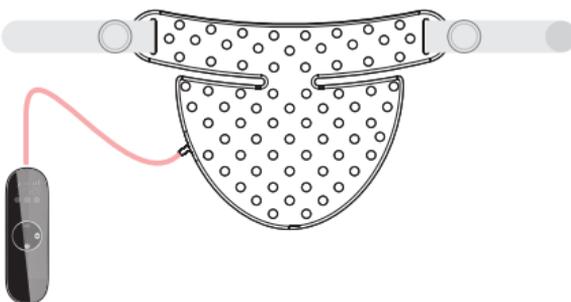Connect the remote control to the mask using the USB-C cable. Connect the second cable to the DC port on the remote for charging.

MODES

Modes	Recommended Serum	Function*
NIR Near-Infrared	Botanical Lifting Effect Serum, Slow-Aging Ampulles, EMS & Light Therapy Serum	Near-Infrared (NIR) light penetrates deeper into the skin layers at 850nm than visible light, providing a warming effect. This can promote cell renewal, improve circulation and skin moisture, and reduce fine lines, wrinkles, and inflammation.
Red Light	Botanical Lifting Effect Serum, Hyaluron Serum, Squalane Serum, EMS & Light Therapy Serum	At 660nm, the Red Light penetrates the dermis, promoting collagen production and vitality, improving skin elasticity and hydration, and imparting a fresh glow.
Yellow Light	Brightening Serum, Vitamin C Serum, EMS & Light Therapy Serum	At 590nm, the Yellow Light accelerates blood circulation, activates cells, balances skin structure, promotes cell growth, and breaks down skin melanin, thereby reducing dark spots, redness, and blue marks.
Blue Light	Bakuchiol Serum, Niacinamide Serum, EMS & Light Therapy Serum	The Blue Light at 410nm is renowned for its antibacterial properties, reducing inflammation and blemishes, and increasing skin oxygenation.

*The benefits described are based on experience.

NOTICE: BEFORE USE

- Before each use, inspect the LED Décolleté Mask. Do not continue use if there are obvious damages.
- The device cannot be operated while charging.
- Store the LED Décolleté Mask flat at room temperature, in dry conditions, and away from direct sunlight.
- Avoid prolonged use of the tool. Excessive use may cause skin issues. It is recommended to use the device 3-5 times per week for 10 to 20 minutes each session. Blue Light mode should not be used for longer than 10 minutes.
- Temporary redness may occur immediately after treatment, which should subside within 1-2 hours.
- If you experience discomfort during use, discontinue immediately and seek medical assistance.
- Keep the LED Décolleté Mask out of the reach of children (0-8 years), who should not perform any maintenance or cleaning on the device. The device is not a toy.
- To protect your eyes, do not stare directly into the LED lights during treatment. Ensure no part of the LED lights cover your eyes during use.
- This product is an electronic device for home use, not a medical device, and is not intended for diagnostic, therapeutic, or commercial purposes.
- Please use only the provided cable with DC port to charge the LED Décolleté Mask. The controller's indicator light will blink during charging.

APPLICATION

BEFORE FIRST USE

Perform a sensitivity test on your décolleté before the first application:

1. Place the mask on your inner arm.
2. Press and hold the ON/OFF button until the lights illuminate, and allow the device to complete its full 10-minute treatment program.
3. If your skin feels hot or uncomfortable at any time, turn off the device immediately and remove it.
4. Wait at least 6 hours after treatment. If redness or itching occurs on your inner arm, the LED Décolleté Mask may not be suitable for you. Consult your doctor before use.

POWER ON

Press and hold the ON/OFF button for approximately 3 sec until the lights illuminate.

SELECTING MODE

Use the ON/OFF button to switch between the three modes by pressing it briefly. The top of the remote control displays an indicator with R | B | Y, showing which mode is active.

Use the NIR button to activate the Near-Infrared mode. This mode can be selected in addition to the three main modes for more intensive treatment.

ADJUSTING INTENSITY

Use the INTENSITY button to adjust the strength. The top of the remote control displays an indicator with three bars indicating the intensity level.

Start with intensity level 1 during your first use to allow your skin time to adjust and build tolerance. Monitor how your skin responds and make adjustments as needed for optimal results.

APPLICATION

SETTING TIME

The treatment starts with a default time of 10 minutes. Use the TIME button to extend the treatment duration. The top of the remote control displays an indicator showing the selected treatment duration of 10, 15, or 20 minutes.

Adjust the treatment duration according to your skin type. Monitor how your skin responds and make adjustments as needed for optimal results.

IMPORTANT: The Blue Light mode should not be used for more than 10 minutes and is preset accordingly.

HOW TO USE

1. Cleanse your décolleté with a gentle cleansing product.
2. Use the mask directly on cleansed skin. Or, if you prefer, apply a serum suitable for light therapy to your skin and allow it to fully absorb.
3. Place the LED Décolleté Mask on your décolleté, loop the straps around your neck, and secure them at the back of your neck with the Velcro. Ensure a comfortable fit.
4. Press the ON/OFF button on the remote control until the lights illuminate to turn on the device. This starts the treatment directly in the Red Light mode, preset to 10 minutes.
5. Switch between modes as desired by briefly pressing the ON/OFF button. Additionally, you can activate the Near-Infrared function with the NIR button.
6. Extend the treatment time as desired by pressing the TIME button.
7. The device will automatically shut off after the set time.
8. Remove the LED Décolleté Mask and proceed with your usual skincare routine. After the treatment, active ingredients penetrate even deeper into the skin.

APPLICATION

For optimal results, use the LED Décolleté Mask 3-5 times per week. Each session should last between 10 and 20 minutes. Following these recommended guidelines maximizes the cosmetic benefits of light therapy for your skin. skin and guarantees a safe usage. The Blue Light mode should not be used for more than 10 minutes.

POWER OFF

The device will automatically turn off after your treatment, depending on the selected time of 10 to 20 minutes. You can end the treatment at any time by pressing the ON/OFF button for 3 seconds.

CLEANING

Clean the LED Décolleté Mask before and after each use. Wipe the mask with a damp cloth or a water-based wipe. Do not submerge the mask in water. Ensure the device is turned off and not charging during cleaning.

BATTERY

The LED Mask operates on a built-in rechargeable lithium-polymer battery, which can be properly recycled at the end of its lifecycle. Always follow correct charging instructions listed below:

1. To charge the mask, attach the supplied DC cable to the remote control and connect it to a power source. The remote control can remain connected to the LED Décolleté Mask. **IMPORTANT:** The mask cannot be charged via the USB-C port.
2. The controller's indicator light will blink during charging. A full charge takes approximately 1 hour.
3. Disconnect the device from the charger after 1 hour of full battery charge to avoid overcharging, which may reduce battery lifespan over time.
4. The device cannot be operated while charging.

BATTERY

5. Only charge the device at ambient temperatures between 0°C/32°F and 40°C/104°F. Avoid charging outdoors, in bathrooms, or within 3 meters of water sources.

CAUTION!

1. The following persons should not use this product:
 - Children (0-8 years)
 - Pregnant or breastfeeding women
 - People with allergies
 - People with skin injuries, eczema, swelling, skin diseases, atopic dermatitis or sensitive skin
 - People with any pre-existing inflammation
 - People with mental illnesses and impaired mental abilities
 - People with a pacemaker or other devices
 - People with fever and other physical abnormalities, heart disease, serious illnesses, acute illnesses, malignant tumours, tuberculosis, blood disorders and other symptoms
 - People with epilepsy or a history of skin cancer (consultation with a doctor is required before use)
 - People with any existing inflammation or pre-existing conditions
2. Do not use on eyes, eyeballs, wounds, plastic parts, anterior neck or larynx, tumours, eczema, over the heart, etc.
3. Do not use the tool in case of any pre-existing inflammations or other illnesses.
4. Avoid using the device after recent injuries, surgery or aesthetic treatments, including neurotoxin treatments, dermal fillers, microneedling, laser and chemical peels, until the skin has fully healed.
5. After aesthetic treatments such as fillers, hyaluronic acid, Botox etc., consult your doctor before use.

CAUTION!

6. If you feel any discomfort during use, discontinue use immediately.
7. The tool should not be used over a long period of time and should always be used in accordance with the instructions for use. Excessive use can cause skin problems. A maximum of 3-5 uses per week for 10-20 minutes each time should not be exceeded. The Blue Light mode should not be used for longer than 10 minutes.
8. The different modes of the LED décolleté mask can generally be used in succession. However, make sure that the maximum total application time of 3-5 per week of 10-20 minutes each, and a maximum of 10 minutes for the Blue Light mode, is not exceeded. If you have sensitive skin, it is also advisable to use the modes on different days to avoid irritation.
9. If you have sensitive skin, consult your dermatologist before first use.
10. Individuals allergic to light or experiencing light allergies from medication should not use this device.
11. If redness or swelling of the skin persists for more than 2 hours, discontinue use and seek medical attention.
12. Do not use the device on open wounds.
13. Do not use the device on painful, elevated, crusted, or bleeding skin. This may be a sign of a more serious condition and should be checked by your doctor.
14. Carry out a sensitisation test on your face as described in the operating instructions before your first application.
15. Use the mask only for the conditions specified in the application guidelines. The device has been tested only for these listed conditions, and risks for other conditions are unknown.
16. Do not apply Retinol before using the Red Light and NIR modes. However, using Bakuchiol is allowed.
17. Some individuals may experience redness, dryness, or irritation after using Blue Light.
18. Avoid excessive exposure of the skin to Blue Light, as it may lead to pigmentation or darkening.

CAUTION!

19. Consult a dermatologist or healthcare provider before using Blue Light if you have an existing skin condition or have undergone dermatological treatments.
20. Your skin may become more sensitive to sunlight for two or three days - apply sunscreen.
21. Do not store the tool near a fire source, heat it, or use it at very high ambient temperatures. Store it in a normally temperate, dry, and ventilated place.
22. If the tool is damaged by impacts or falls, please do not continue to use it.
23. Do not modify, disassemble, or repair the tool yourself. For repairs, please contact the distributor or manufacturer directly.
24. Prevent the internal battery from coming into contact with skin, eyes, or mouth. Wear protective gloves if exposed to battery fluids and seek medical help if skin or eye contact occurs.
25. Do not use accessories not included in the package.
26. Do not twist, bend, forcibly pull out the charging cable, wrap it around the device, or otherwise damage it, as this may cause fire or short circuit. If the charging cable is damaged or becomes hot, stop using it immediately.
27. Avoid operating the device under blankets or pillows to prevent excessive heating and potential fire hazards.
28. Do not operate the device without removing creams, serums, lotions, essential oils, and makeup before treatment. These products may cause photosensitivity or affect the effectiveness of the LED treatment. Apply serums suitable for light therapy before treatment.

CARE AND MAINTENANCE

1. The LED Décolleté Mask and accessories are not waterproof and should not be submerged in water or machine-washed. Avoid placing the device where it can fall into or be drawn into water.

CARE AND MAINTENANCE

2. Stains and product residues on the mask can be wiped off with a damp cloth after use.
3. Do not use alkaline cleaners, volatile oils with alcohol, thinners, or other solvents for cleaning.
4. If the tool is not used for an extended period, it should be powered on regularly and charged every three months to check functionality. Clean the tool before longer periods of non-use and store it in a sealed, dry place.
5. Do not open the device, batteries, or charger, as this poses risks of electric shock, fire, or exposure to battery chemicals. Disassembly or improper reassembly will void the warranty.
6. Use only the specified cleaning products to avoid damage.
7. Store the LED Décolleté Mask, battery, and charging cable in a cool, dry place.

DISPOSAL

IMPORTANT: Return the tool to a authorized collection point for disposal after battery is completely empty. Do not disassemble the tool yourself.

This symbol indicates that batteries and accumulators must not be disposed of with regular household waste. You are legally obligated to return them free of charge to a publicly authorized collection point or to collection points set up by distributors for recycling at the end of their lifespan.

Batteries and accumulators that are not firmly enclosed in the device and can be removed without destruction must be separated and disposed of in a discharged state according to the designated disposal method. The recycling and proper disposal of old devices and batteries are important contributions to protecting our environment. Improper disposal can lead to toxic substances

DISPOSAL

entering the environment, which can have harmful effects on humans, animals, and plants.

Below are some best practices and recommendations to extend the lifespan of your batteries and accumulators and enable their reuse:

- Charge your batteries properly and fully to maximize their lifespan. Fully discharge them with a suitable charger if necessary before recharging.
- Consider using rechargeable batteries. These can be reused and recharged multiple times, reducing the amount of battery waste.
- Always use the correct type of batteries for your devices. Improper usage can shorten battery life and may have harmful effects.

Always follow safety instructions when handling used batteries. Lithium-containing used batteries (Li = Lithium) pose a high fire hazard. Therefore, special attention must be paid to the proper disposal of lithium-containing used batteries and accumulators. Improper disposal can also lead to internal and external short circuits due to thermal effects (heat) or mechanical damage. A short circuit can result in a fire or explosion and have serious consequences for humans and the environment. To prevent a short circuit, tape over the terminals or contacts of the battery before disposal.

WARRANTY

1. The product is covered under warranty for proper use as per instructions. Parts and repairs will be offered free of charge within one year from the date of purchase for malfunctions or damages.
2. Product defects or damages caused by human error, misuse, or improper use are not covered by the warranty.

WARRANTY

3. Strong collisions or impacts, excessively high or low voltage power supply, as well as product defects or damages caused by unavoidable circumstances, are not covered by the warranty.
4. For the above cases not covered by the warranty, the distributor will decide if repair is possible and charge for parts and repairs accordingly.
5. This warranty clause applies only to products purchased and used by consumers for personal use, not for industrial and commercial use.
6. We reserve the right to adjust the product image due to product optimization. The product corresponds to the actual object.
7. Manufacturer & Distributor: Rosental Organics GmbH, Rosa-Luxemburg-Str. 14, 10178 Berlin, Germany

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