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
## USER MANUAL

---

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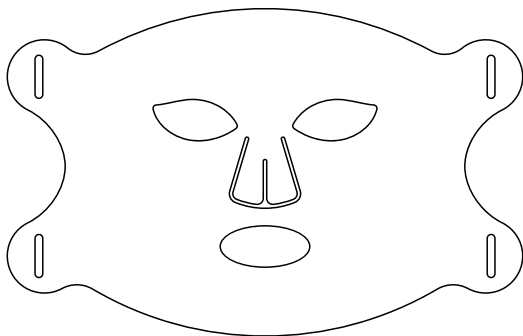**We wish you lots of enjoyment with your  
LED Face Mask | Advanced Beauty Tool  
with light therapy from Rosental Organics!**

Before using your new tool, please read this manual carefully and keep it for future reference.


To ensure correct usage of this product, please read all safety precautions carefully before using. These precautions will help you use this product correctly to avoid personal injury or damage to yourself or others. The safety precautions are listed below. Please comply with the corresponding regulations!


# PARTS

---

## Remote Control

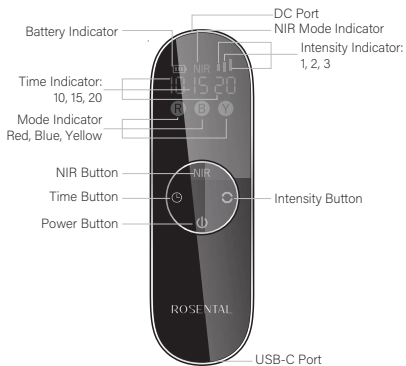
## LED Face Mask
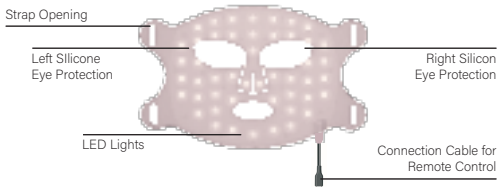
# ASSEMBLY

---

## STEP 1

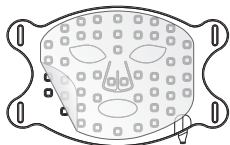Remove the protective film from the mask.

## STEP 2

Attach the right and left silicone eye protection inside the mask.

---


## STEP 3

Attach the straps to the mask by threading them through the openings and securing them inside with the Velcro.

---

## STEP 4

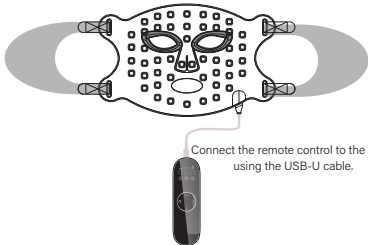Connect the remote control to the mask using the USB-U cable.

## MODES

---

Mode	Recommended Serum	Function
<b>NIR Near-Infrared</b>	Skin Activation Serum, Botanical Lifting Effect Serum, Slow-Aging Ampullen	Near-Infrared (NIR) light penetrates deeper into the skin layers at 850nm than visible light, providing a warming effect. This can promote cell renewal, improve circulation and skin moisture, and reduce fine lines, wrinkles, and inflammation.
<b>Red Light</b>	Skin Activation Serum, Botanical Lifting Effect Serum, Hyaluron Serum, Squalane Serum	At 660nm, the Red Light penetrates the dermis, can promote collagen production and vitality, improving skin elasticity and hydration, and imparting a fresh glow.
<b>Yellow Light</b>	Skin Activation Serum, Brightening Serum, Vitamin C Serum	At 590nm, the Yellow Light can accelerate blood circulation, activates cells, balances skin structure, promotes cell growth, and breaks down skin melanin, thereby reducing dark spots, redness, and blue marks.
<b>Blue Light</b>	Skin Activation Serum, Bakuchiol Serum, Niacinamide Serum	The Blue Light at 410nm can have anti-bacterial properties, reducing inflammation and blemishes, and increasing skin oxygenation.

## NOTICE: BEFORE USE

---

- Before each use, inspect the LED Face Mask. Do not continue use if there are obvious damages.
- The device cannot be operated while charging.
- Store the LED Face Mask flat at room temperature, in dry conditions, and away from direct sunlight.
- Avoid prolonged use of the tool. Excessive use may cause skin issues. It is recommended to use the device 3-5 times per week for 10 to 20 minutes each session.
- Temporary redness may occur immediately after treatment, which should subside within 1-2 hours.
- If you experience discomfort after use, discontinue immediately and seek medical assistance.
- Keep the LED Face Mask out of reach of children. This device is not a toy.
- To protect your eyes, do not stare directly into the LED lights and use eye protection during treatment. Ensure no part of the LED lights cover your eyes during use.
- This product is an electronic device for home use, not a medical device, and is not intended for diagnostic, therapeutic, or commercial purposes.
- Please use only the provided cable to charge the LED Face Mask. The controller's indicator light will blink during charging.

## APPLICATION

---

### BEFORE FIRST USE

Perform a sensitivity test on your face before the first application:

1. Place the mask on your inner arm.
2. Press and hold the ON/OFF button until the lights illuminate, and allow the device to complete its full 10-minute treatment program.
3. If your skin feels hot or uncomfortable at any time, turn off the device immediately and remove it.
4. Wait at least 6 hours after treatment. If redness or itching occurs on your inner arm, the LED Face Mask may not be suitable for you. Consult your doctor before use.

### POWER ON

Press and hold the ON/OFF button for approx. 3 sec. until the lights illuminate.

### SELECTING MODE

Use the ON/OFF button to switch between the three modes by pressing it briefly. The top of the remote control displays an indicator with R | B | Y, showing which mode is active.

Use the NIR button to activate the Near-Infrared mode. This mode can be selected in addition to the three main modes for more intensive treatment.

### ADJUSTING INTENSITY

Use the INTENSITY button to adjust the strength. The top of the remote control displays an indicator with three bars indicating the intensity level. Start with intensity level 1 during your first use to allow your skin time to adjust and build tolerance. Monitor how your skin responds and make adjustments as needed for optimal results.

## APPLICATION

---

### SETTING TIME

The treatment starts with a default time of 10 minutes. Use the TIME button to extend the treatment duration. The top of the remote control displays an indicator showing the selected treatment duration of 10, 15, or 20 minutes. Adjust the treatment duration according to your skin type. Monitor how your skin responds and make adjustments as needed for optimal results. The Blue Light mode should not be used for more than 10 minutes and is preset accordingly.

### HOW TO USE

1. Cleanse your face with a gentle facial cleanser.
2. Apply a serum suitable for light therapy to your skin.
3. Place the LED Face Mask on your face, loop the straps around your ears, and secure them at the back of your head with the Velcro. Ensure a comfortable fit. Press the ON/OFF button on the remote control until the lights illuminate to turn on the device. This starts the treatment directly in the Red Light mode, preset to 10 minutes.
4. Switch between modes as desired by briefly pressing the ON/OFF button. Additionally, you can activate the Near-Infrared function with the NIR button.
5. Extend the treatment time as desired by pressing the TIME button.
6. The device will automatically shut off after the set time.
7. Remove the LED Face Mask and proceed with your skincare routine.

For optimal results, use the LED Face Mask 3-5 times per week. Each session should last between 10 and 20 minutes. Following these recommended guidelines maximizes the benefits of light therapy for your skin.

## APPLICATION

---

### POWER OFF

The device will automatically turn off after your treatment, depending on the selected time of 10 to 20 minutes. You can end the treatment at any time by pressing the ON/OFF button for 3 seconds.

### CLEANING

Clean the LED Face Mask before and after each use. Wipe the mask with a damp cloth or a water-based wipe. Do not submerge the mask in water. Ensure the device is turned off and not charging during cleaning. You can remove the eye attachments and clean them under running water. Dry them with a cloth afterward.

## BATTERY

---

The LED Mask operates on a built-in rechargeable lithium-polymer battery, which can be properly recycled at the end of its lifecycle. Always follow correct charging instructions listed below:

1. To charge the mask, attach the provided DC cable to the remote control and connect it to a power source. The remote control can remain connected to the LED Face Mask during charging.
2. The controller's indicator light will blink during charging. A full charge takes approximately 1 hour.
3. Disconnect the device from the charger after 1 hour of full battery charge to avoid overcharging, which may reduce battery lifespan.
4. The device cannot be operated while charging.
5. Only charge the device at ambient temperatures between 0°C/32°F and 40°C/104°F. Avoid charging outdoors, in bathrooms, or within 3 m of water sources.

## CAUTION!

---

1. The following individuals should not use this product:
  - Children, pregnant or breastfeeding women
  - Individuals with allergies
  - People with skin injuries, eczema, swelling, skin diseases, atopic dermatitis, or sensitive skin
  - Individuals with mental illnesses
  - People with pacemakers or other devices
  - Individuals with fever and other physical abnormalities, heart diseases, severe illnesses, acute diseases, malignant tumors, tuberculosis, blood disorders, and other symptoms
2. Do not use on eyes, eyeballs, wounds, plastic parts, front neck or larynx, tumors, eczema, above the heart, etc.
3. Avoid using the device after recent injuries, surgeries, or facial treatments, including neurotoxin treatments, dermal fillers, microneedling, laser and chemical peels, until the skin is fully healed.
4. If you have sensitive skin, consult your dermatologist before first use.
5. Individuals allergic to light or experiencing light allergies from medication should not use this device.
6. If redness or swelling of the skin persists for more than 2 hours, discontinue use and seek medical attention.
7. Do not use the device on open wounds.
8. Do not use the device on painful, elevated, crusted, or bleeding skin. This may be a sign of a more serious condition and should be checked by your doctor.
9. Use the mask only for the conditions specified in the application guidelines. The device has been tested only for these listed conditions, and risks for other conditions are unknown.
10. Do not apply Retinol before using the Red Light and NIR modes. However, using Bakuchiol is allowed.

## CAUTION!

---

11. Some individuals may experience redness, dryness, or irritation after using Blue Light.
12. Avoid excessive exposure of the skin to Blue Light, as it may lead to pigmentation or darkening.
13. Consult a dermatologist or healthcare provider before using Blue Light if you have an existing skin condition or have undergone dermatological treatments.
14. Your skin may become more sensitive to sunlight for two or three days - apply sunscreen.
15. Do not store the tool near a fire source, heat it, or use it at very high ambient temperatures. Store it in a normally temperate, dry, and ventilated place.
16. If the tool is damaged by impacts or falls, please do not use it.
17. Do not modify, disassemble, or repair the tool yourself. For repairs, please contact the distributor or manufacturer directly.
18. Prevent the internal battery from coming into contact with skin, eyes, or mouth. Wear protective gloves if exposed to battery fluids and seek medical help if skin or eye contact occurs.
19. Do not use accessories not included in the package.
20. Do not twist, bend, forcibly pull out the charging cable, wrap it around the device, or otherwise damage it, as this may cause fire or short circuit. If the charging cable is damaged or becomes hot, stop using it.
21. Avoid operating the device under blankets or pillows to prevent excessive heating and potential fire hazards.
22. Do not operate the device without removing creams, serums, lotions, essential oils, and makeup before treatment. These products may cause photosensitivity or affect the effectiveness of the LED treatment. Apply serums suitable for light therapy before treatment.

## CARE AND MAINTENANCE

---

1. The LED Face Mask and accessories are not waterproof and should not be submerged in water or machine-washed. Avoid placing the device where it can fall into or be drawn into water.
2. Stains and product residues on the mask can be wiped off with a damp cloth after use.
3. Do not use alkaline cleaners, volatile oils with alcohol, thinners, or other solvents for cleaning.
4. If the tool is not used for an extended period, it should be powered on regularly and charged every three months to check functionality. Clean the tool before longer periods of non-use and store it in a sealed, dry place.
5. Do not open the device, batteries, or charger, as this poses risks of electric shock, fire, or exposure to battery chemicals. Disassembly or improper reassembly will void the warranty.
6. Use only the specified cleaning products to avoid damage.
7. Store the LED Face Mask, battery, and charging cable in a cool, dry place.

## DISPOSAL

---

1. Before disposal, dismantle the tool properly and remove the battery from the remote control. First, turn off the power supply.
2. Carefully disconnect the battery connector from the circuit board and remove the battery.
3. Dispose of the battery separately from the housing properly at a recycling station according to local guidelines.

## WARRANTY

---

1. The product is covered under warranty for proper use as per instructions. Parts and repairs will be offered free of charge within one year from the date of purchase for malfunctions or damages.
2. Product defects or damages caused by human error, misuse, or improper use are not covered by the warranty.
3. Strong collisions or impacts, excessively high or low voltage power supply, as well as product defects or damages caused by unavoidable circumstances, are not covered by the warranty. For the above cases not covered by the warranty, the distributor will decide if repair is possible and charge for parts and repairs accordingly.
4. This warranty clause applies only to products purchased and used by consumers for personal use, not for industrial and commercial use.
5. We reserve the right to adjust the product image due to product optimization. The product corresponds to the actual object.
6. Manufacturer & Distributor: Rosental Organics GmbH,  
Rosa-Luxemburg-Str. 14, 10178 Berlin, Germany

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