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USER MANUAL

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
USER MANUAL

Thank you for purchasing Rosental Organics' EMS Gua Sha | Beauty Tool with Light Therapy!

Before using your new Face Yoga Tool, please read this manual carefully and keep it for future reference.

To ensure that you use this product correctly, please read all safety precautions carefully before use. The safety precautions can help you use this product correctly to avoid personal injury or property damage to you or others. The safety precautions are as follows, please follow the relevant regulations of the safety precautions!

COMPONENTS

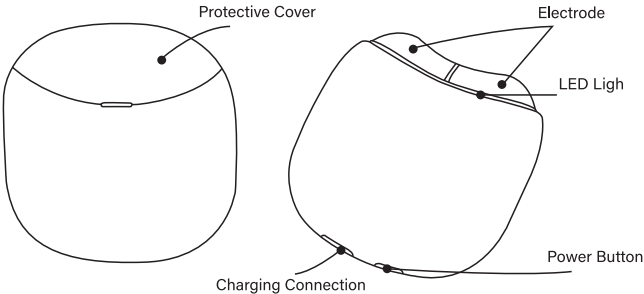BLUE LIGHT	GREEN LIGHT	RED LIGHT
Clarifying	Balancing	Slow-Aging
Low EMS	Middle EMS	Strong EMS
▪ anti-inflammatory ▪ reduces the formation of bacteria ▪ counteracts impurities	▪ relieves redness and tension ▪ promotes a smooth teint ▪ balancing	▪ supports collagen production ▪ provides more firmness and elasticity ▪ stimulates blood circulation and metabolism

NOTICE: BEFORE USE

- Before each use of this tool, please check whether there is any obvious damage. If it is damaged, please do not use it.
- The tool should not be used continuously for a long time. Excessive use may cause skin problems. A usage of maximum 5 minutes each time, 2-3 times a week is recommended.
- The skin's sensitivity to heat compress and EMS tolerance are different. Users with very sensitive skin may experience redness or swelling after prolonged use. Please adjust the usage time according to individual needs.
- If you feel any discomfort after use, please stop using it immediately and seek medical assistance.
- Please keep it out of reach of children.
- Excessive use of this product may damage the skin.
- Do not look directly at the LED light to avoid damage to the eyes.
- Due to the EMS function of the product, an electric shock like feeling on the skin is normal.
- This product is a home use electronic product, not a medical device, and is not intended for diagnostic, therapeutic and other commercial purposes.

HOW TO USE

POWER ON

Press and hold the button for about 2 seconds, the product is turned on, and the LED light, EMS and warming function start to run. For first-time use, please be sure to use the Blue Light mode with lowest EMS strength to make sure your skin can adapt to the device.

FOREHEAD DESTRESSING

Sweep from the center of the forehead to the ears.

DEFINED CHEEBONES

Sweep from the upper and lower lip up to the ear. Then from the tip of the chin jawbone to the ear.

JAWLINE AND NECK TIGHTENING

Sweep from the lower jaw across the cheeks to the upper jaw. The same from the bottom of the jaw along the neck to the collarbone.

CHANGE EMS MODE

To switch the EMS mode shortly press the power button.

POWER OFF

To switch off your tool, press and hold the power button for about 2 seconds to shut down. The tool will automatically shut down after 5 minutes of continuous use in the same mode.

CLEANING

Before and after each use, please wipe the electrode of the gua sha with a clean soft cloth dipped in a little water to clean it.

BATTERY INSTRUCTION

This tool is powered by a built-in rechargeable lithium polymer battery, which can be properly recycled at the end of its lifecycle. Please always follow the correct instructions for charging the battery, as listed below:

1. When the battery of the product is low, there will be a low voltage alarm: the red LED light will flash and automatically shut down.
2. When the battery voltage is low, connect with the charging cable and connect the USB plug at the other end to a power outlet. Please use the original charging cable to charge this tool.
3. The red LED light flashes while the product is charging, the green LED light is on when fully charged. The green LED light will turn off after 30 minutes if the charging cable is not disconnected after being fully charged. It takes about 1.5 hours to fully charge.

ATTENTION!

1. The following people should not use this tool:
 - children
 - pregnant or menstruating or breastfeeding women
 - People with allergies
 - People with skin trauma, eczema, swelling, skin disease, atopic dermatitis or sensitive skin
 - People with mental illness
 - People with a cardiac rhythm regulator or other device
 - Fever and other physical abnormalities, heart disease, serious illness, acute illness, malignant tumor, tuberculosis, blood disease and other symptoms

2. Do not use on the eyes, eyeballs, wounds, plastic parts, neck or laryngeal nucleus, above the heart, etc.
3. Do not immerse this product in water or other liquids.
4. Do not use or place the tool close to a fire source, heat it, or use it in a high temperature environment. Please store the tool in a normal temperature, dry and ventilated place.
5. If the product is damaged due to bumping or falling, please stop usage immediately.
6. Do not modify, disassemble, or repair the tool yourself. For repairs, please contact the distributor or manufacturer directly.
7. Keep the tool away from strong electromagnetic fields, such as induction.
8. Do not use accessories that are not included in the scope of delivery.
9. Do not twist, bend, forcefully pull out, wrap the charging cable around the device, or otherwise damage it, as this may cause fire or short circuit. If the charging cable is damaged or becomes hot, stop using it immediately.
10. Do not use with essential oils. Exclusively use with serums or oil developed for face yoga.

CARE AND MAINTENANCE

1. The tool is not waterproof and should not be washed with or immersed in water.
2. Stains and product residue on the electrode head and product can be wiped off with a damp cloth after use.
3. Do not use alkaline cleaning agents, volatile oils containing alcohol, thinners, or other solvents for cleaning.
4. When the tool is not in use for an extended period, it should be regularly powered on and charged every three months to check its functionality. Before prolonged non-use, clean the tool and store it in a sealed and dry environment.

DISPOSAL

1. Before disposing of the tool, properly disassemble it and remove the battery. First, turn off the power of the tool.
2. Use a flat screwdriver to carefully loosen the bottom casing of the tool.
3. Gently disconnect the battery plug from the circuit board and remove the battery.
4. Dispose of the battery separately from the casing at a designated recycling station.

TROUBLESHOOTING

If a malfunction occurs during operation, please check the following. If the fault cannot be rectified or the product still cannot operate ideally, please contact the customer support directly via contact@rosental.de and do not disassemble and attempt to repair it by yourself.

FAULT PHENOMENOM	CAUSE ANALYSIS	TROUBLESHOOTING
The device can not be turned on.	The battery voltage is low, and the self-protection function is activated.	Please fully charge the device before use.
The device can not be charged.	The charge cable or adapter is not connected to power.	Connect the charging cable correctly and make sure the adapter is on.
The device is powered off after red LED flashes.	The battery voltage is low.	Please use the device after fully being charged.

WARRANTY

1. The product is in normal use according to the instructions, if it is faulty or damaged, parts and maintenance will be free of charge within one year from the date of purchase.
2. Product failure or damage caused by human error and abuse or misuse is not covered by the warranty.
3. Strong collisions or clashes, too high or too low power supply voltage, and product failure or damage caused by irresistible reasons are not covered by the warranty.
4. For the above-mentioned situation not covered by warranty, the distributor will determine whether it can be repaired and collect parts and maintenance fees.
5. This warranty clause is only applicable to occasions purchased and used by consumers for personal consumption needs, not for industrial and commercial use.
6. We keep the right to amend the product picture due to product upgrading. The product is subject to the real object.
Manufacturer: Shenzhen Margotan Tech Co., Ltd
7. Address: 12F, Building #A, Rongchao Binhai Mansion, No.2021 Haixiu Road, Bao'an Central District, Shenzhen, China.
8. Distributor: Rosental Organics GmbH

SPEZIFIKATION / SPECIFICATION

Produktname | Product Name

EMS Gua Sha | Beauty Tool with Light Therapy

Produktmodell | Product Model

IF-GS02

Ladespannung/-strom | Charging Voltage/Current

5 V / 0,5 A

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Polymer-Lithium-Batterie 3,7V 600mAh | Polymer lithium battery 3.7V 600mAh

Nennleistung | Rated Power

≤ 4 Watt

Produktmaterial | Product Material

ABS

Produktabmessungen | Product Measurements

93,2 x 82 x 27,6 mm

Produktgewicht | Product Weight

90 g

Herkunftsland | Country of Origin

PRC

ROSENTHAL
ORGANICS